


MMT ACCESSORIES FOR UROFLOWMETRY ARE EASY TO USE AND MEET THE HIGH STANDARDS OF DEMANDING HOSPITAL ENVIROMENTS



Reusable stainless steel container sterilisable in an autoclave



Plastic containers with or without measuring scale

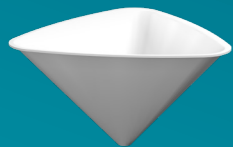

Reusable funnel, 121°C resistant



Disposable paper funnel

Notify body No. 2265, 3EC International a.s.

MMTD_2

MEDKONSULT medical technology s.r.o.
Pasturova 15, 779 00 Olomouc, Czech Republic

T +420 581 113 030
E mmt@mmtsystems.com

MMT DANFLOWS ARE COMPATIBLE WITH A WIDE RANGE OF ACCESSORIES

Three different types of the stand designed to fit underneath the commode chair for seated voiding.



Standard stand

Small stand designed for videourodynamics

Portable version with a folding stand



Uroflowmetry stand is easy to open/lock and adjust the height.



Lightweight and easy to clean commode chair also suitable for videourodynamic procedure.



GET INTO THE FLOW WITH MMT URODYNAMICS



Danflow

www.mmtsystems.com

www.mmtsystems.com

YOUR CHOICE ON REPORTING AND ARCHIVING UROFLOWMETRY EXAMINATIONS



UROFLOWMETRY DEVICE CONNECTED TO A COMPUTER OR NOTEBOOK VIA CABLE OR WIRELESSLY.

Advanced SW functionality includes automatic/manual uroflowmetry, patient database, real time flow curves, multiple types of nomograms.

SMALL COMPACT THERMAL PRINTER CONNECTED WITH THE FLOW METER VIA CABLE OR WIRELESSLY. Automatic data storage onto an inbuilt SD card. SD memory card can be easily removed and inserted into your computer for subsequent storage or processing records if desired.

FLOW METERS WITH EITHER WIRED OR WIRELESS WEIGHT CELL TRANSDUCER TECHNOLOGY TO CREATE CLASSIC UROFLOWMETRY RECORDINGS/ANALYSIS



Cable connection



Unique technology at 868/915 MHz

Massively improved range between flow rate and printer unit or PC



Wireless technology Bluetooth



The waterproof flow meter allows easy cleaning and disinfection. The lithium battery provides incredible battery performance for a truly wireless experience. No need for daily charging.

WITH THE DANFLOW RANGE OF URINE FLOW METERS, MMT HAVE MADE UROFLOWMETRY SIMPLE AGAIN!



BASIC FEATURES OF FLOW METERS DANFLOW

EASY TO OPERATE

EASY TO INSTALL

METHODS CONFORM TO ICS STANDARDS

Basic features

	DANFLOW Chord CU	DANFLOW Wave CU	DANFLOW Chord PC	DANFLOW Wave
Nomograms	-	-	•	•
Real time curves	-	-	•	•
Automatic / manual mode	•	•	•	•
Patients' database	-	-	•	•
Post-processing	-	-	•	•
PDF protocol	•	•	•	•
Customised protocols	-	-	•	•
Language versions	•	•	•	•
Online remote service support	-	-	•	•
Dicom	-	-	•	•
HL7	-	-	•	•
Wired connection	•	-	•	-
Wireless connection - 868/915 MHz	-	•	-	•
Wireless connection - Bluetooth	-	•	-	•
Thermoprinter	•	•	-	-
UFM stand - standard	•	•	•	•
UFM stand - small	•	•	•	•
UFM stand - portable	•	•	•	•
Commode chair	•	•	•	•
Reusable / disposable accessories	•	•	•	•

• available - not available