

COST EFFECTIVE SYSTEM COVERING ALL STANDARD URODYNAMIC EXAMINATIONS



UNIQUE WIRELESS CONNECTION
AT 868/915 MHZ ENSURES STABILITY
AND NO SIGNAL DROPOUT

COMPATIBILITY WITH DIFFERENT TYPES
OF THE PRESSURE MEASURING SYSTEMS

Notify body No. 2265, 3EC International a.s.

MMTM_2

MEDKONSULT medical technology s.r.o.
Pasteurova 15, 779 00 Olomouc, Czech Republic

T +420 581 113 030
E mmt@mmtsystems.com



UROMIC Melody is a small compact system for complete urodynamic examinations.

The unit is equipped with 5 measuring channels, an output for a profilometric puller and inbuilt cystometric pump.

UROMIC Melody enables wireless or USB connection with a laptop or standard PC.

Available accessories meet the high standards of the demanding clinical environment and enable easy handling during the examination.



UROMIC Melody unit can be mounted either on a table stand or on an infusion stand and easily moved into position and convenient for compact storage.

The infusion stand is specially designed for use with MMT urodynamics systems and includes a second pole to allow full range transducer adjustment as well as infusion bag hooks.

Accessories can be fixed at any height and easily removed.

Braked wheels ensure patient and system safety.

www.mmtsystems.com


GET INTO THE FLOW WITH MMT URODYNAMICS



UROMIC *Melody*

www.mmtsystems.com

USER FRIENDLY SOFTWARE UDMvision

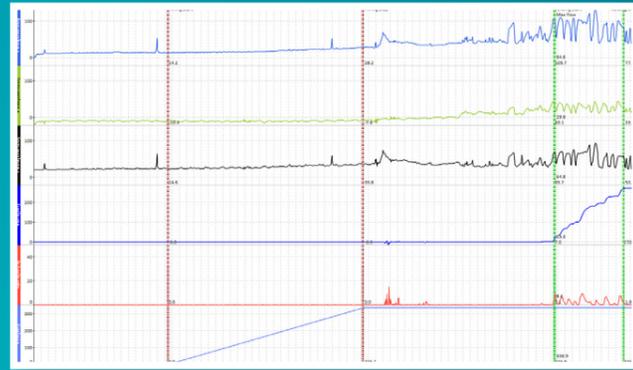

EASY TO INSTALL

PREDEFINED PROTOCOLS

METHODS CONFORM TO ICS STANDARDS

MULTIPLE NOMOGRAMS

CONFIGURABLE EXAMINATION METHODS AND TEST REPORTS

ANORECTAL MANOMETRY

VIDEOURODYNAMICS

PVR VALUE AUTOMATICALLY SAVED INTO TEST REPORTS

CONFIGURABLE LANGUAGES

WIRELESS CONNECTION WITH HOSPITAL NETWORK, INCLUDING DICOM AND HL7

VOICE OR REMOTE CONTROL

ONLINE REMOTE SERVICE SUPPORT



MMT URODYNAMIC SYSTEMS ARE COMPATIBLE WITH A WIDE RANGE OF ACCESSORIES



Melady

WIRELESS PATIENT UNIT



Portable, wireless UROMIC Harmony provides additional wireless channels for use at the patient side, keeping cables out of the way.

PROFILOMETRY



MMT unique magnetic puller with customisable speed and resterilisable catheter holder.

VALSALVA LEAK POINT PRESSURE*



Flexible camera allows easy visualisation of urethral leaks.

*Third party product

CYSTOMETRY



MMT cystometric pump and infusion weight ensure precise volume filling.

EMG / BIOFEEDBACK



EMG curve is smooth and precise thanks to an external preamplifier.

3D BLADDER SCANNER*



The portable, intuitive, and accurate bladder scanner provides real-time PVR information. The PVR value is available in seconds and saved into test reports automatically.

*Third party product

UROFLOWMETRY



The waterproof flow meter allows easy cleaning and disinfection. The lithium battery provides incredible battery performance for a truly wireless experience. No need for daily charging.



Different types of stand for flow meters available - height adjustable, designed to fit underneath the commode chair. Easy to assemble folding version allows easy travelling and more flexibility.



Lightweight and easy to clean commode chair for seated voiding, suitable also for videourodynamics procedures.